

Schedule DAY 2

TUESDAY	ACTIVITY	PRESENTER
9.00 – 9.10	Review of Day 1	Jean
9.10 – 10.00	Lab: Amplicon QC	Louise
10.00 – 10.30	Morning tea	
10.30 – 12.30	Lab: Library Preparation & Loading	Jean
12.30 – 13.30	Lunch	
13.30 – 14.15	Lab: Flow cell loading practice and Run Reviewing	Louise
14.15 – 14.30	Lecture: REVIEW OF LABORATORY WORKFLOW	Jean
14.30 – 15.30	Panel: Ask us anything!	
15.30 – 16.00	Afternoon Tea	
16.00 – 17:00	Lecture: MPXV Illumina Sequencing Options	Jean

REVIEW OF LABORATORY WORKFLOW

Jean Moselen
Senior Medical Scientist
Victorian Infectious Disease Laboratory (VIDRL)

Genomics key goals

—

Performing NGS experiments accurately and reproducibly.

Requires a careful and considerate mode of working.

NGS is a highly sensitive technique.

=

Highly sensitive to **contamination**

=


Incorrect patient diagnosis or **invalid** run.

Mpox surveillance workflow.

The core stages of this workflow are:

- a) specimen collection
- b) sample preparation
- c) genome sequencing processing of sequencing results
- d) sequence data interpretation and data sharing.

Location of work in the NGS lab is important

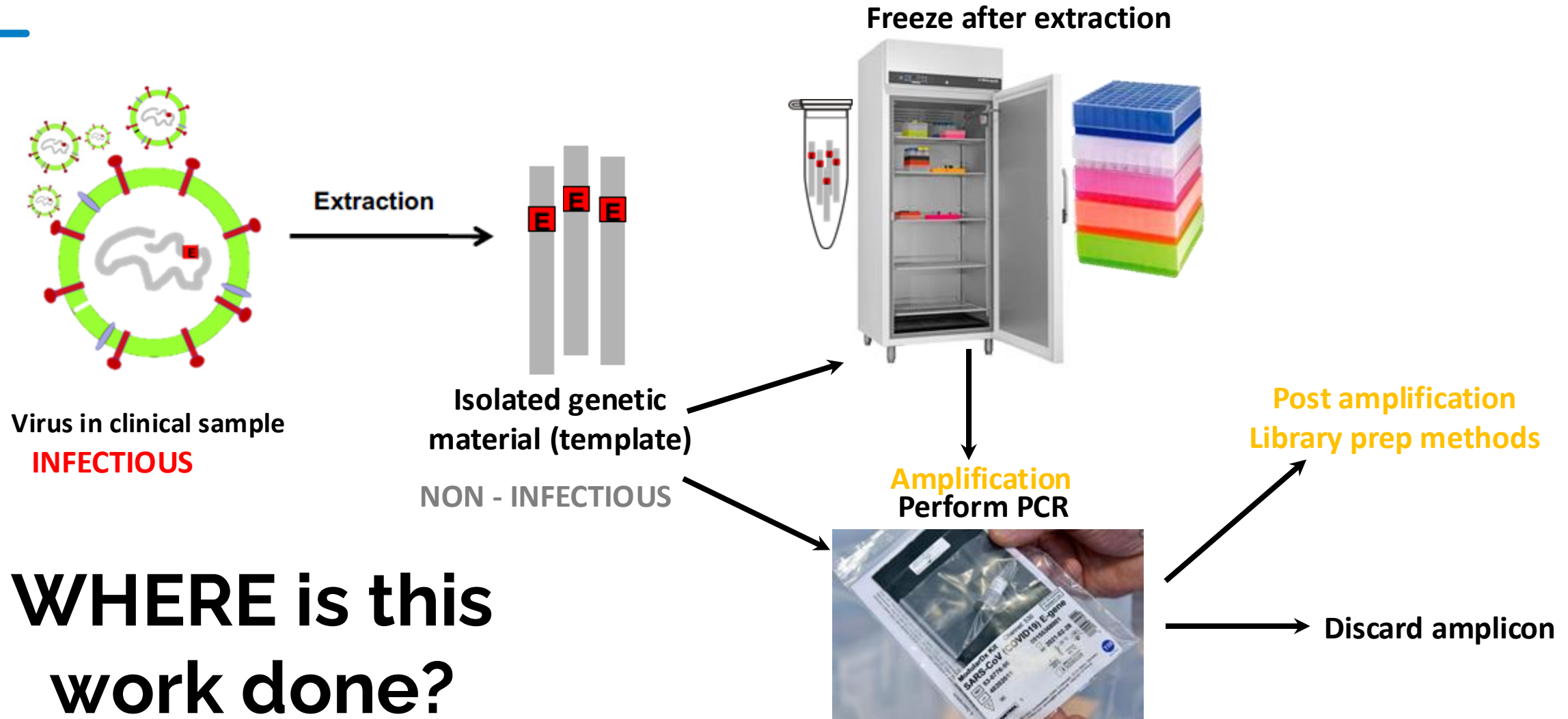


Next Generation Sequencing

Zones of work

Work can be done over days/weeks

Safe stopping points

Emphasis is always on the **avoidance of contamination** with each step of the workflow...

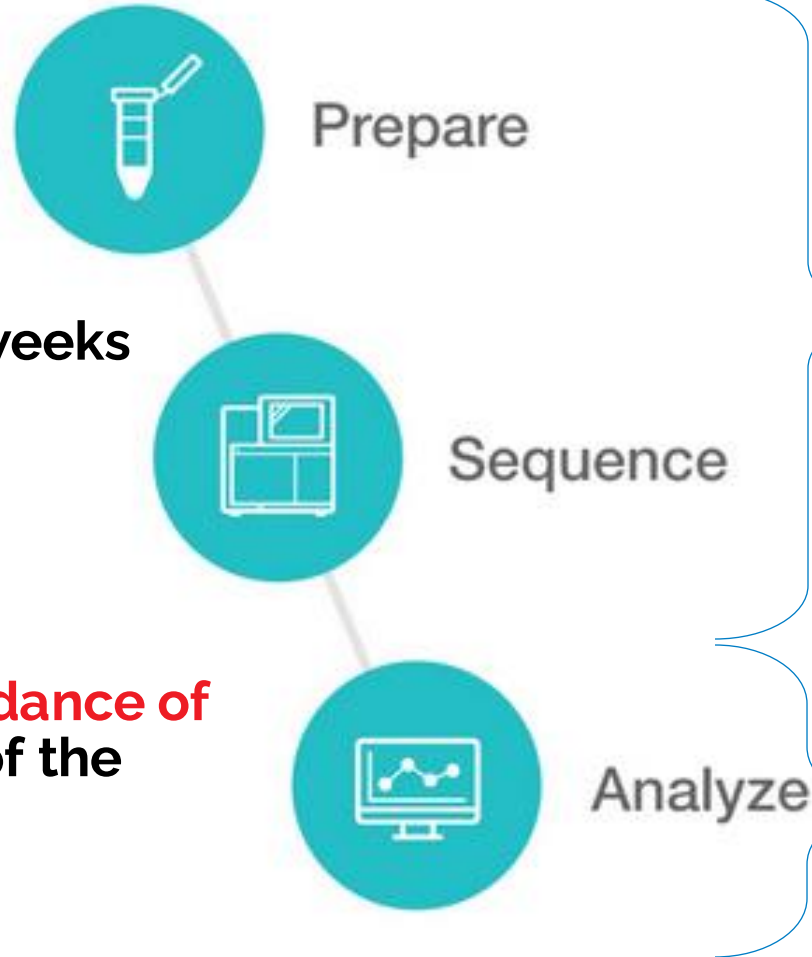


WET LAB

DRY LAB

Any kind of NGS technology generates a significant amount of data.

UNIDIRECTIONAL WORKFLOW

ONE DIRECTION 1 → 2 → 3 → 4 → 5 → 6

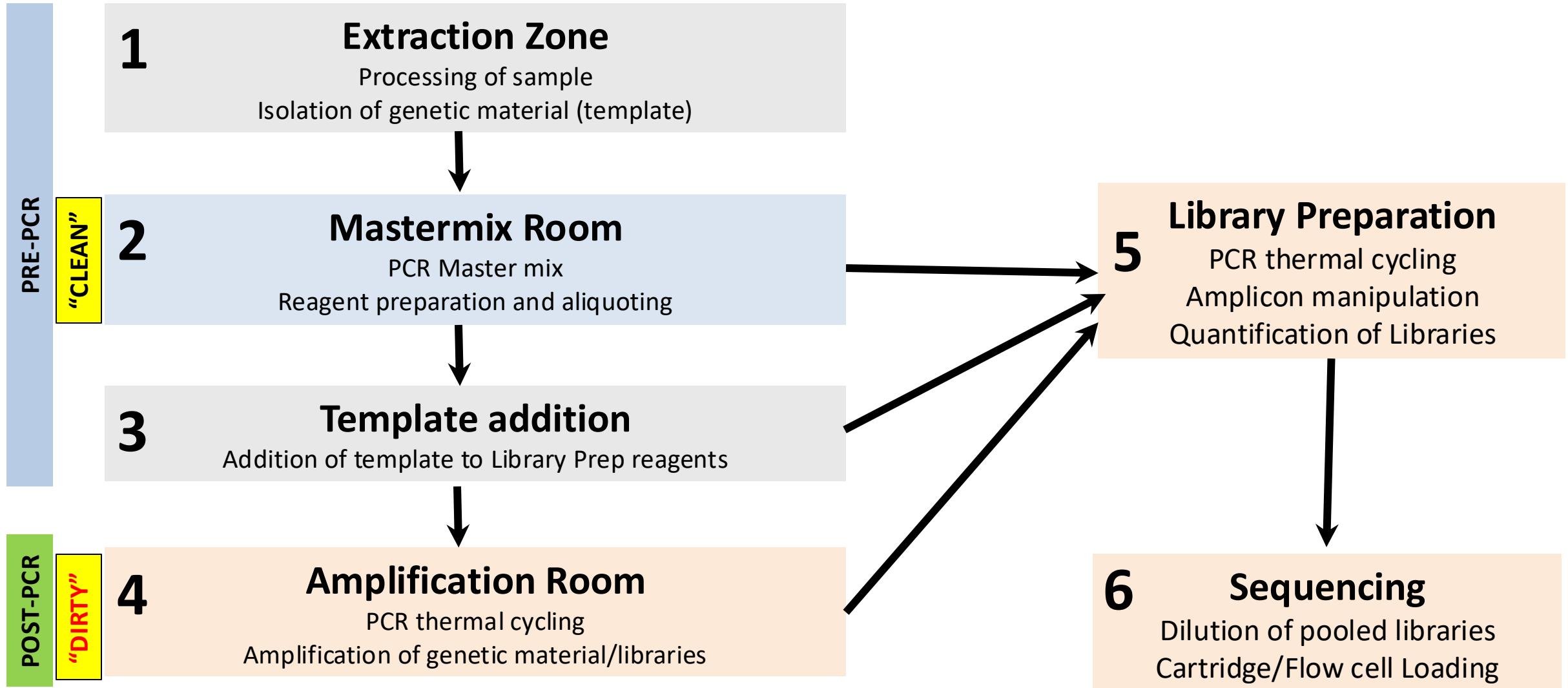

NGS lab Infrastructure

Wet lab – Equipment vs throughput

HIGH <96 samples +

96 well plate format

EXTRACTION



King Fisher Apex
+Tan Bead Maelstrom™ 9610

QUANTIFICATION




CLARIOstar Plus

ANALYSER




TapeStation 4200

LIQUID HANDLER

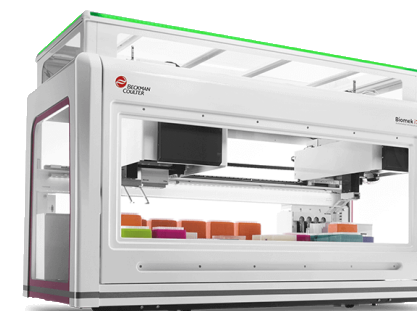

Biomek i5 x 3 & i7 x 2

LOW <48 samples

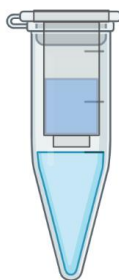


Column extractions



QuBit 4



TapeStation 2200

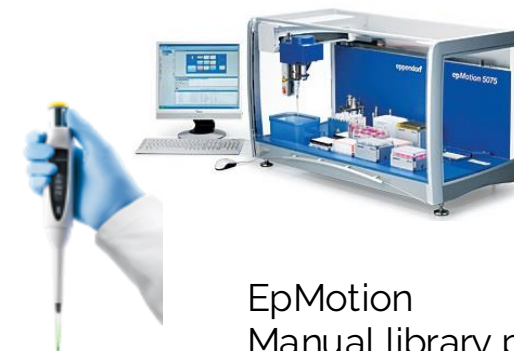

EpMotion
Manual library prep

NGS lab Infrastructure

Wet lab – Sequencer vs throughput

Data Output low(er)

GridION



iSeq100



Amplicon schemes + Bait capture low sample #

MiSeq



Amplicon schemes
Bait capture
Complex samples

Data Output high(er)

NextSeq 1000/2000



Deep Sequencing
= lots of data

PromethION 2 Solo



Metagenomics
Bait capture
Complex samples
Adaptive sampling

Routine cleaning

Equipment and surfaces

All areas during lab procedures (particularly lab bench, PCR workstation and BSCII)

Scheduled routine cleaning. Daily/weekly basis.

Use a DNAAWAY followed by 80% ethanol to wipe down your:

Bench top – Spray/wipe

Pipettes – wipe only

Centrifuge and Vortex – wipe only

Racks – Spray/soak (DET-SOL 5000)

Freezer Doors, cabinets, drawers – spray/wipe

Thermocycler lid and buttons - wipe

Genomics contamination


The best way to deal with the risk of contamination is to take every precaution to prevent it from happening!

Environmental checks

contamination by exogenous DNA or RNA must be controlled to minimize the risk of faulty results.

AMPLIFIED MATERIAL CARRYOVER

Products from previous NGS
Not following unidirectional workflow

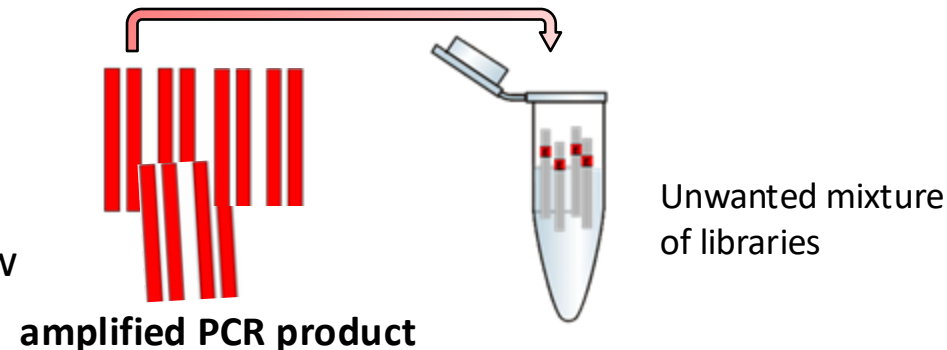


Contamination in Library Prep

Opening more than one tube at a time
OR
Not changing gloves after a spill or pipetting
error
OR
Not changing tips between samples

can lead to....

Sample prep area,
extraction, during library
prep (pre barcoding)

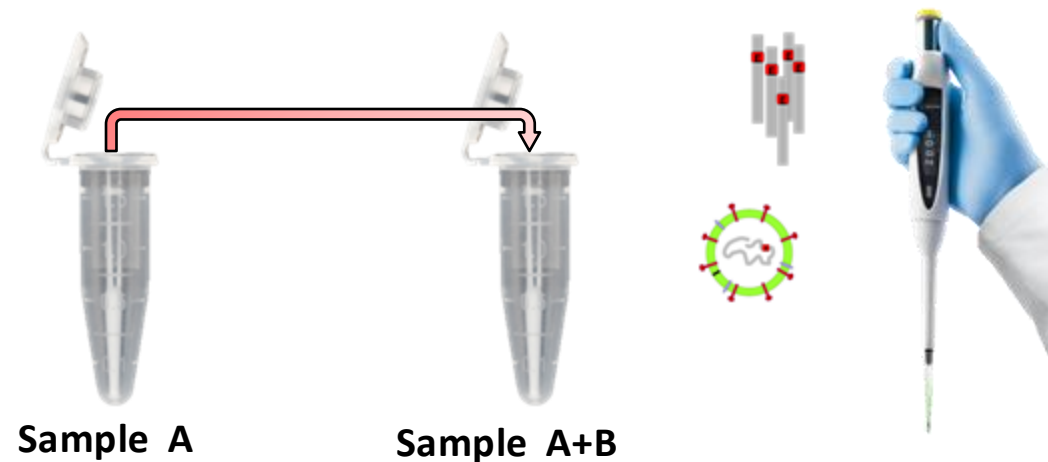

Storage of NGS samples and reagents

Samples & Extracted material / cDNA

- -20°C (freeze/thaw considerations)
- “to be sequenced” box

Completed Library (plate) or
Pooled Library (tube)

- Freeze if not used immediately

Data

- Have backup hard Drives, S Drives etc.



Sequencing with ONT

What can you sequence on ONT?

Strand sequencing method = nanopore sequences what ever is presented to it regardless of length

Any nucleic acid – Genomic = bacterial, fungal, metagenome

PCR amplicons – Tiled SARS CoV- 2, RSV, universal amplicons for segmented viral genomes e.g. flu.

PCR panels - AMR detection

RNA – Viral genomes, transcriptomes, cDNA

Targeted sequencing = cas9 enrichment, hybridization enrichment, adaptive sampling.

Sample Input

Library Preparation kits have input requirements

- Sample volume, concentration of genetic material, sample type

Accurate quantification of nucleic acid is critical for NGS methods



Fluorometric (Qubit)– Targets specific type nucleic acid using a fluorophore and detected by a fluorometer, greater sensitivity. Less impacted by contaminants Qubit.
Binds = specific detection of molecules e.g. RNA or DNA

Sample pretreatment

Library Preparation kits have input requirements

GOAL= greater proportion of pathogen reads in library


Pathogens can make up a small proportion of the original specimen

- DNase treatment
- Ampure bead clean up (1:1)
- Ribosomal RNA removal
- **Nucleic Acid concentration** (Kits, elution volume, speed vac)
- **cDNA Synthesis** (First strand and Second Strand)



highly concentrated, purified ssDNA/RNA that is suitable
for subsequent library prep methods

Bioinformatics

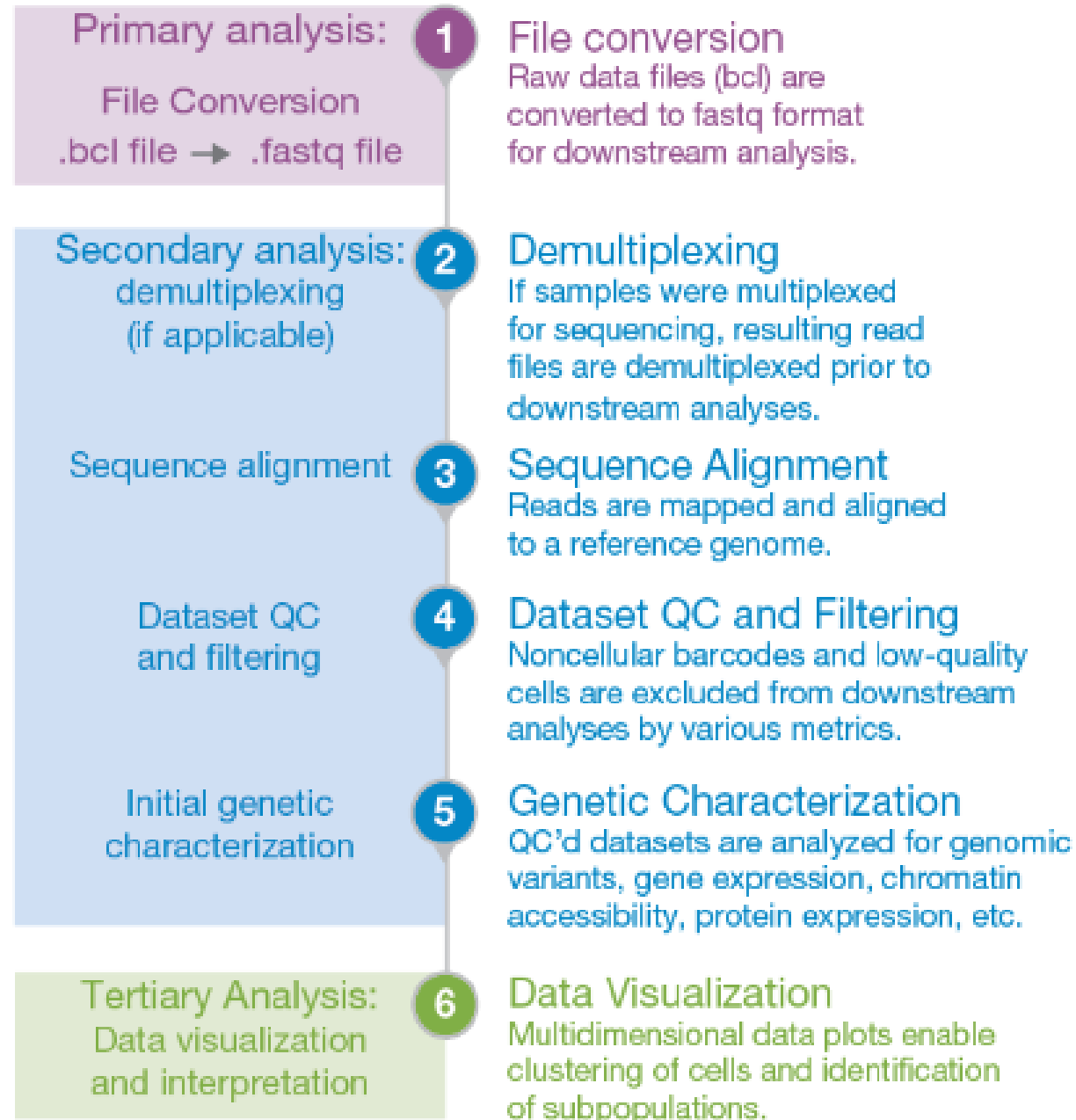

Schedule DAY 2

TUESDAY	ACTIVITY	PRESENTER
9.00 – 9.10	Review of Day 1	Jean
9.10 – 10.00	Lab: Amplicon QC	Louise
10.00 – 10.30	Morning tea	
10.30 – 12.30	Lab: Library Preparation & Loading	Jean
12.30 – 13.30	Lunch	
13.30 – 14.15	Lab: Flow cell loading practice and Run Reviewing	Louise
14.15 – 14.30	Lecture: Review of Laboratory Workflow	Jean
14.30 – 15.30	Panel: Ask us anything!	
15.30 – 16.00	Afternoon Tea	
16.00 – 17:00	Lecture: MPXV Illumina Sequencing Options	Jean