

# Wet Lab schedule DAY 1

MONDAY	ACTIVITY	PRESENTER
8:45 – 9.00	Registration	Lisa
9.00 – 9.15	Overview of the Doherty Institute/CPG/VIDRL/MDU	Lisa
9.15 – 9:45	Welcome and introductions	
9:45 –10.00	Training overview	Jean
10.00– 10.30	LAB: Pipetting exercise	Louise
10.30 – 11.00	Morning tea	
11.00 – 11.30	LECTURE: Introduction to MPXV genomics and VIDRL	Jean
<b>11.30 – 12.30</b>	<b>LECTURE: TILED AMPLICON FOR MPOX</b>	<b>Jean</b>
12.30 – 13.30	Lunch	
13.30 – 15.30	LAB: Tiled amplicon PCR	Louise
15.30 – 16:00	Afternoon tea	
16:00 – 16.30	LECTURE: Introduction to ONT sequencing viruses	Louise
16.30 – 17:00	Group discussion: Opportunity for Q&A and further discussion	Nicole

# WET LAB DAY 1: TILED AMPLICON

---

Jean Moselen  
Senior Medical Scientist  
Victorian Infectious Disease Laboratory (VIDRL)

# Viral detection & sequencing methods

## Real time PCR

Basics of using PCR as a **detection method** are very similar for different viruses

## Sanger sequencing

Sequencing only the **gene/s of interest** for identification of pathogen & it's strain/lineage

## Whole-genome sequencing (WGS)

WGS is a comprehensive method for **analysing entire genomes**.

## Next Generation Sequencing (NGS)

**Technology** for sequencing any human, animal, plants, or pathogen microbes.


# Mpox surveillance workflow.

The core stages of this workflow are:


- a) specimen collection
- b) sample preparation – AMPLICON GENERATION**
- c) genome sequencing – ONT PLATFORM**
- d) processing of sequencing results
- e) sequence data interpretation and data sharing.

# VIRAL SEQUENCING METHODS


## Direct metagenomic sequencing

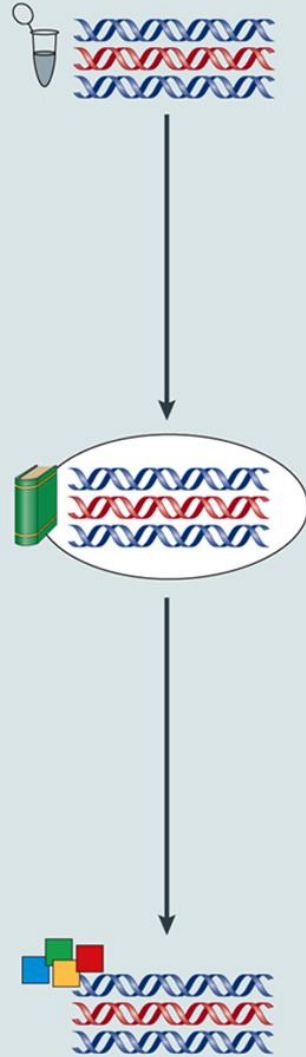

## PCR amplicon sequencing

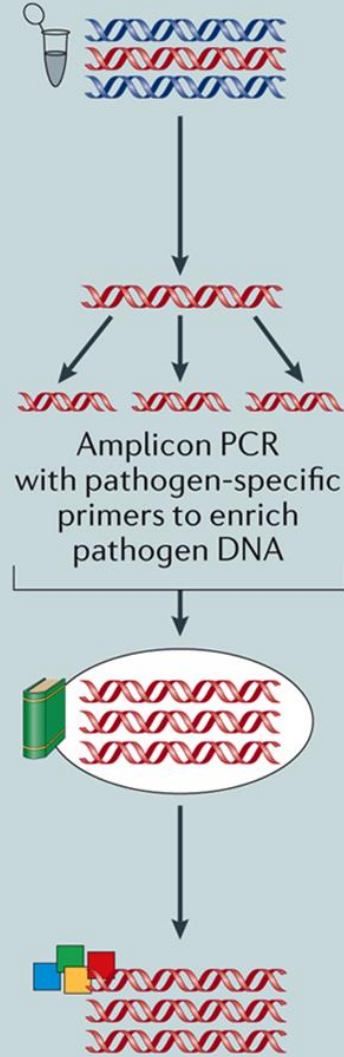


## Target enrichment sequencing

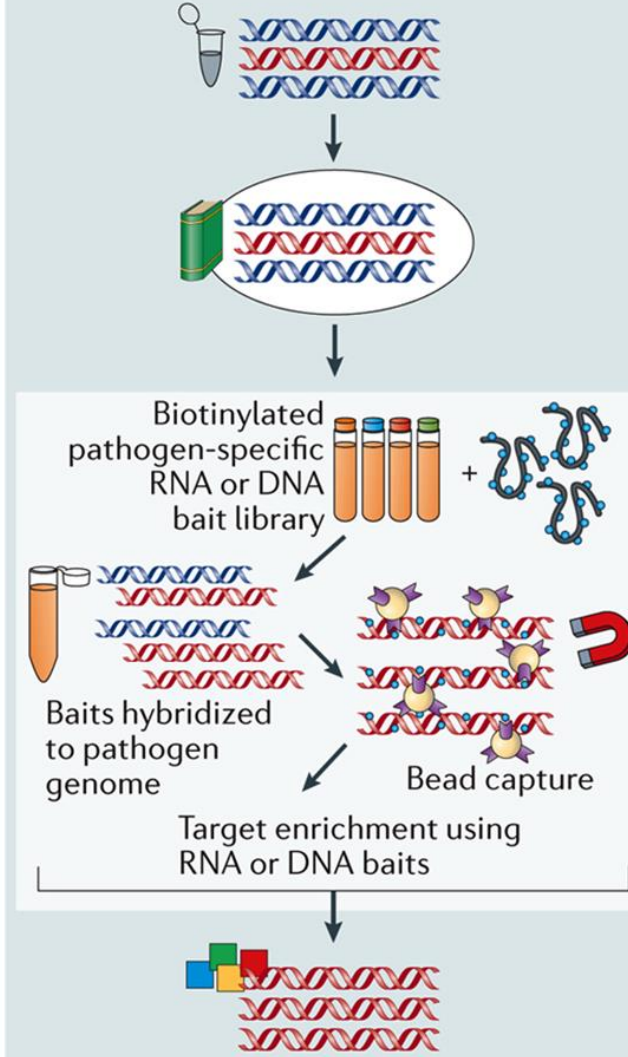

 Pathogen DNA or cDNA

 Host DNA or cDNA

 Clinical sample

 Library preparation

 Sequencing

# Pathogen NGS

## Why so many options?

—

**Metagenomics**

**Amplicon**

**Target Enriched**

Pathogens make up a small proportion of the sample

Viral Diversity

“Needle in the haystack” scenario

Variance in sequencing capacity of laboratories

Different research goals and throughput

Bioinformatic capacity

Emerging field and new technologies

# NGS Schemes Mpox Training

**Abundance of host and other non-target sequences  
in NGS libraries results in low sensitivity**

**How can we get the best result?**

**GOAL= improve # of target specific reads**



Prepare

**Sample type**

**Extraction method**

**Concentration/cleanup**

**Library Input**



Sequence

**Metagenomics**

**Amplicon**

**Target Enriched**

**Sequencer platform**



Analyze

**Bioinformatics tools**


**Reference sequences**

**In house vs online pipelines**




# Viral NGS

## Technique considerations

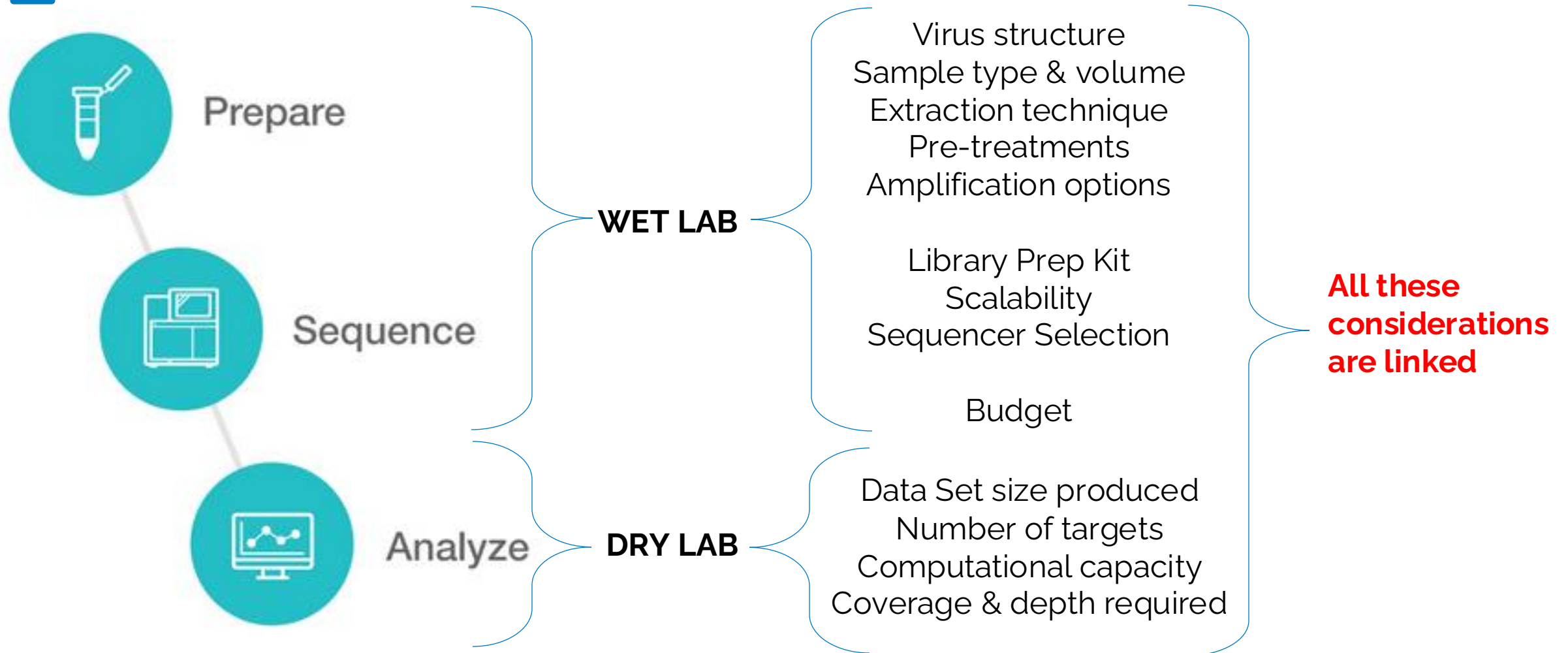

# Viral Genomics Methodology

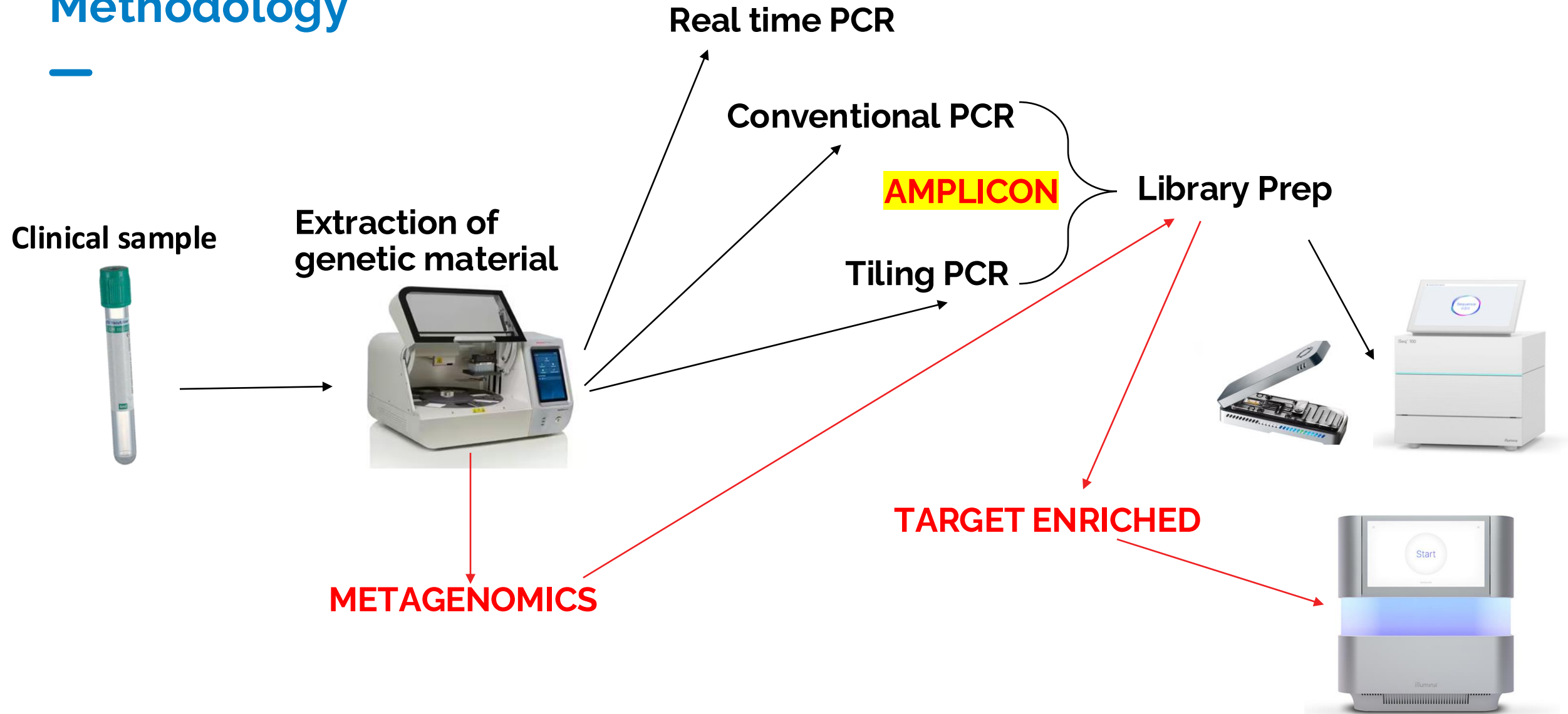

# Amplicon-based WGS

---

Developed to support genomic surveillance efforts.

Used in different public health emergencies, Zika virus and Ebola.

Broadly applied during the SARS-CoV-2 pandemic.

Primary method in viral genomic epidemiology with clinical case-derived samples due to its:

- high sensitivity
- multiplexing capability
- cost-effectiveness
- scalability

# Pan Approach vs Targeted Scheme Design

---

**GOAL= improve # of pathogen specific reads**

**Single target = Greater specificity and sensitivity.**

Higher cost and not useful for detection of dual infections

VS

**PAN approach = multiple targets (viruses) of interest**

e.g Flavi viruses, MMR, respiratory

“PANELS” of viruses tested at once

Useful when looking at among human, animal, and environmental specimens.

Sensitivity is lower

# Outbreak schemes

## Amplicon WGS

---

Workflows were optimized for **simple adaptation of workflows** in laboratories previously involved in SARS-CoV-2 genomic surveillance.

Amplicon-based WGS approaches are highly sensitive due to targeted amplification of the viral DNA, making them particularly suitable to analyze **clinical, environmental and wastewater** samples with low viral loads (up to Ct ~32).

**Labs can order in primer sets and start testing within a short time frame.**

# Sequencing Approaches Mpox

EXAMPLE	SEQUENCING APPROACH
Clinical samples with low viral load (Ct > 28)	Amplicon-based WGS OR Target capture
Detection of structural variants, ITR variation	Amplicon-based WGS or Metagenomic sequencing with long-read technology
Accurate SNP calling and clade typing	Amplicon-based WGS
Validation of amplicon dropout or suspected deletions	Metagenomic sequencing or Target capture
High-throughput routine surveillance	Amplicon-based WGS
Field deployment with minimal infrastructure	Amplicon-based WGS with long-read technology
Wastewater or environmental samples	Metagenomic sequencing or Target capture with short amplicons Possibly: amplicon-based sequencing with short amplicons
Resource-limited labs with basic capacity	Amplicon-based WGS

# Community Designed Amplicon Schemes


 protocols.io



VERSION 2  
FEB 01, 2023

OPEN  ACCESS

DOI:  
[dx.doi.org/10.17504/protocols.io.kqdg39xxeg25/v2](https://dx.doi.org/10.17504/protocols.io.kqdg39xxeg25/v2)

 DengueSeq: A pan-serotype whole genome amplicon sequencing protocol for dengue virus V.2

Chantal Vogels<sup>1</sup>, Chrispin Chaguza<sup>1</sup>, Mallery I Breban<sup>1</sup>,  
Emma Taylor-Abigail J.  
Afeez Sodeinde<sup>1</sup>, Salmon<sup>1</sup>, Porzucek<sup>1</sup>,  
Nathan D Grubaugh<sup>1</sup>

<sup>1</sup>Department of Epidemiology of Microbial Diseases, Yale School of Public Health



Nathan D Grubaugh  
Department of Epidemiology of Microbial Diseases, Yale Schoo...

## ABSTRACT

### Version 2 updates:


- Updated the DENV3 primer file (DENV3\_Primer-Scheme.xlsx) with correct version



FEB 23, 2022

 SHARE

WORKS FOR ME 1

 Poliovirus direct detection and nanopore sequencing (DDNS) FAQs

DOI

[dx.doi.org/10.17504/protocols.io.b5ggq3tw](https://dx.doi.org/10.17504/protocols.io.b5ggq3tw)

Alex Shaw<sup>1</sup>, Catherine Troman<sup>1</sup>, Joyce Akello<sup>1</sup>, Erika Bujaki<sup>2</sup>  
Manasi Majumdar<sup>2</sup>, Javier Martin<sup>2</sup>, Nick Grassly<sup>1</sup>

<sup>1</sup>Imperial College London;

<sup>2</sup>National Institute for Biological Standards and Control

Poliovirus Sequencing Consortium



Nick Grassly  
Imperial College London



# Sequencing Enrichment Approaches for Viruses

## Tiled Amplicon PCR

- Multiple primers designed that span the entire genome --> except 3' and 5' region
- Primer pools to allow overlapping regions amplification while reducing any unwanted primer interactions
- Simple to use
- Lower cost
- Widely used







# Advantages of Tiled Amplicon

- Relatively low cost to sequence
- Works on low viral copy numbers
- Simple workflows
- Faster turn-around-times for results
- Short amplicons (e.g. ~400 bp) more likely to sequence even when genomic material is degraded
- If a region fails to amplify, less coverage is lost compared to if the amplicon were longer
- **Works directly from clinical samples**




# Disadvantages of Tiled Amplicon

- Unable to sequence ends of viral genome --> slightly less complete genome coverage
- Mutations in primer binding sites can cause amplicon drop-outs
- PCR step can amplify small amounts of cross-contamination
- Not ideal for highly diverse or recombinant viruses or new emergence
- Creation of the Primer pools = prone to user error = lots of primers!





## A Whole Segment PCR

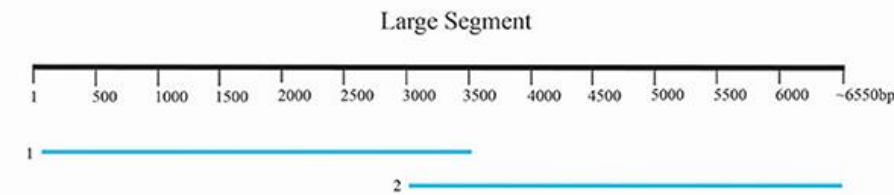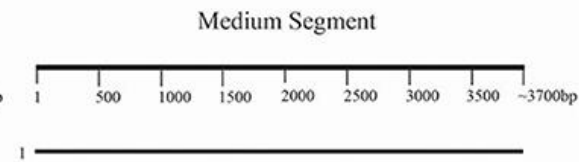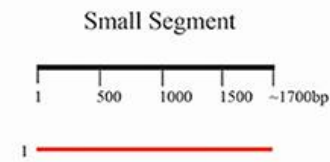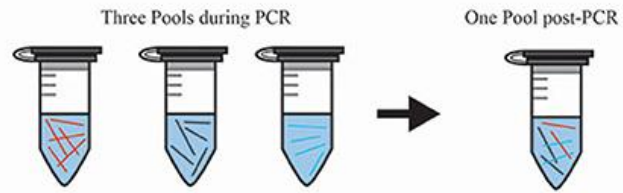





## B Whole Genome Multiplex Tiling

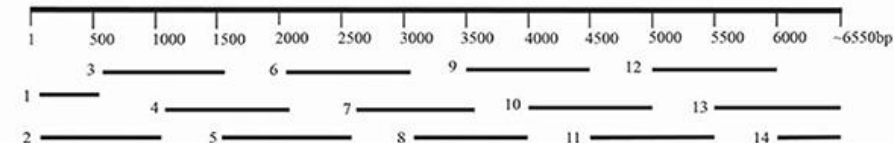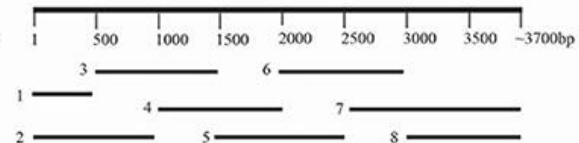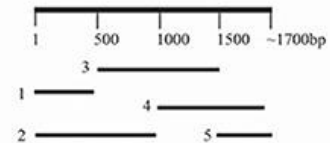

## C Segment-Specific Tiling

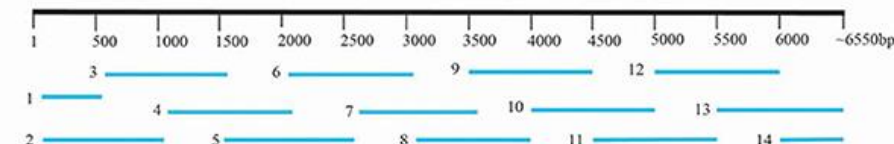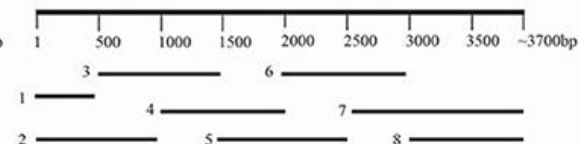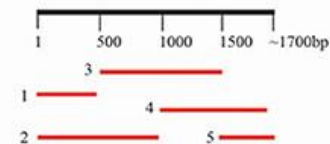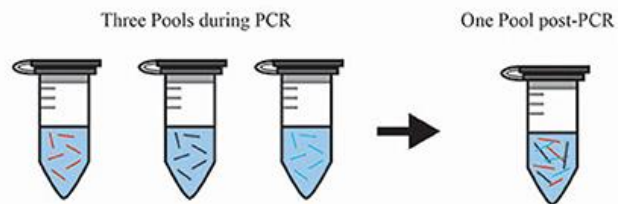

## D Disjointed Tiling

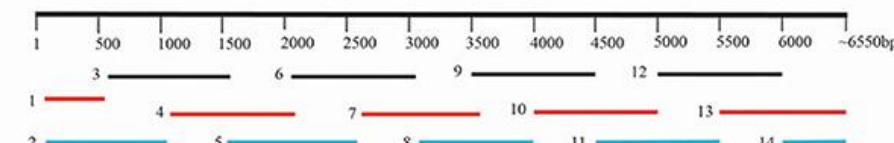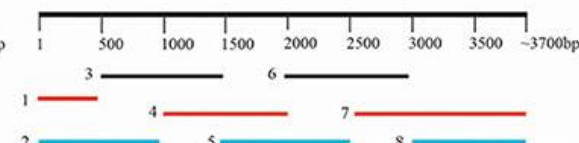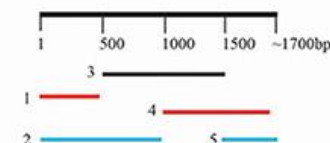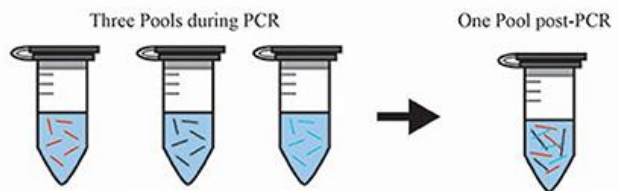

# Sequencing Enrichment Approaches for Viruses

## Tiled Amplicon PCR


Forward and Reverse Primers

Range of dilutions and concentrations

Combinations of primers

Different amplification reactions

= "Pools" of PCR master mixes

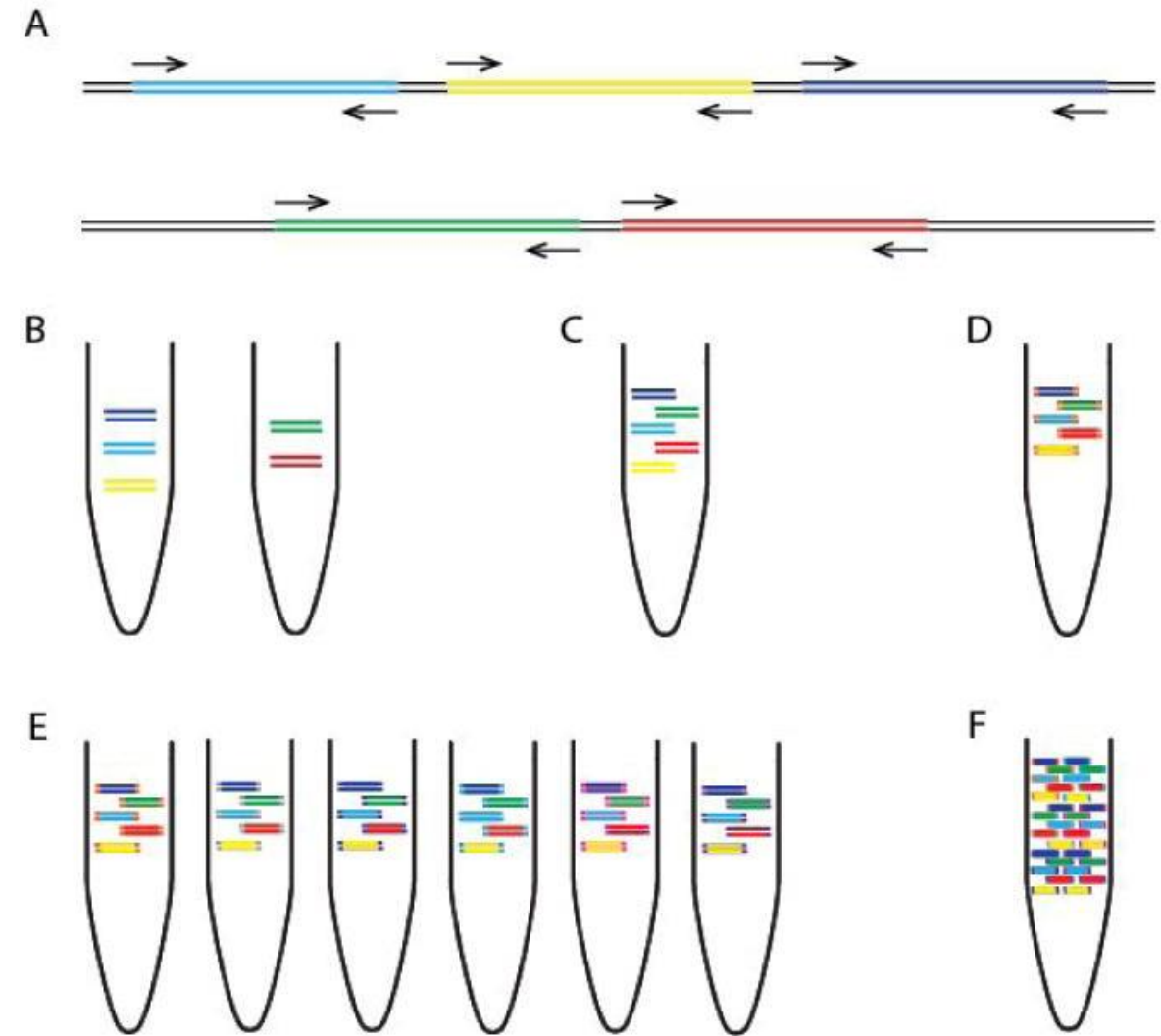

# Mpox Amplicon schemes

## WGS of MPXV

	Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design	Sequencing Platform (validated)	Clade Ib coverage
★	Chen et al. 2023	Average 200	163	<31	IIb (MT903345)	Illumina	85-90% (43)
	Brinkmann et al. 2024	375	682	<30	IIb (ON585033.1)	ONT	
★	Welkers et al. 2022	2500	88	<25	IIb (ON563414.3)	ONT	>93.5 %
	Bosmeny et al. 2023	3000	73	<32	IIb (NC_063383.1)	ONT	
	Isabel et al. 2023	5000	43	<27.9	IIb (ON563414.3)	Illumina	
	ARTIC / INRB	400 (wastewater surveillance) ★ 2500 (clinical surveillance)	147	<30	I and II (KJ642613.1)		
	ARTIC / BCCDC	5000	98		IIb (ON563414.3)		
	Yinda et al. 2023	12,500		<~27		ONT	

★ VIDRL have used these schemes

# Mpox Amplicon schemes

## WGS of MPXV

---

Majority of MPXV primer schemes were developed using reference genomes from Clade IIb resulting in mismatches when applied to Clade I/Ib genomes.

A pan-clade primer scheme has been developed for optimized primer binding across clades

Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design
ARTIC / INRB	400 (wastewater surveillance) 2500 (clinical surveillance)	147	<30	I and II (KJ642613.1)

# Details of Amplicon Schemes

Named via publication  
or authors or lab group

Pool 1 = 76 primers

Pool 2 = 71 primers


Important to know for data  
analysis and bioinformatics

Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design
ARTIC / INRB	400 (wastewater surveillance) 2500 (clinical surveillance)	147	<30	I and II (KJ642613.1)

same scheme name but different  
amplicon options/uses

Guideline for ideal sample quality

# Amplicon Schemes Updates

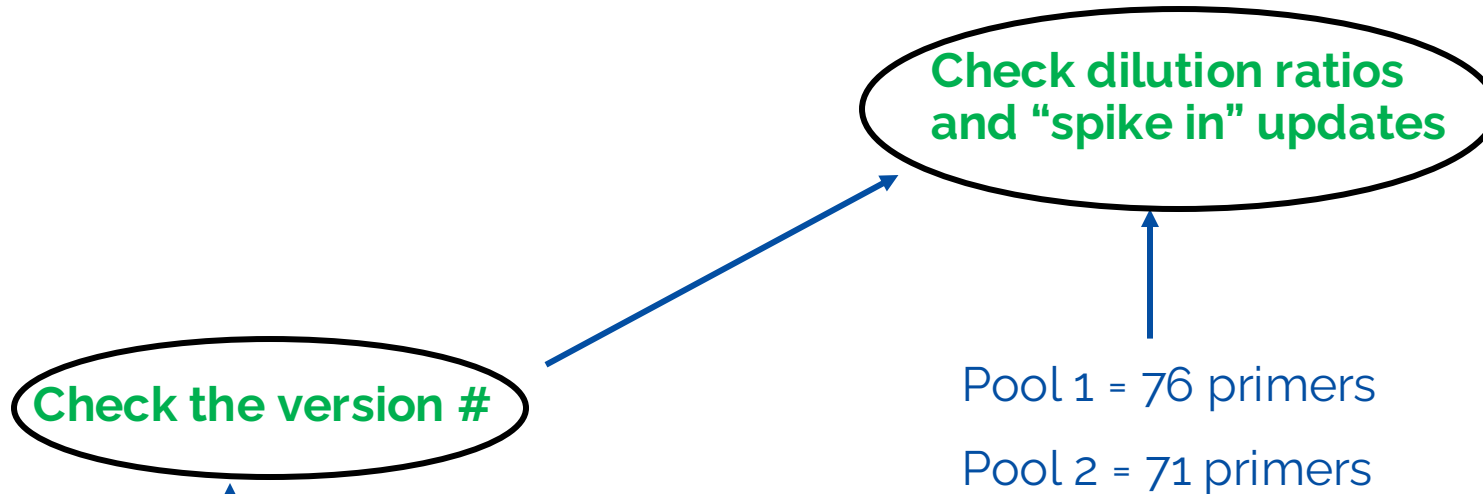

Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design
ARTIC / INRB	400 (wastewater surveillance) 2500 (clinical surveillance)	147	<30	I and II (KJ642613.1)



# Amplicon Schemes Updates



primalscheme labs

[Search](#) [FAQs](#) [About](#) [GitHub](#)

[Create Scheme](#)

## Scheme Details

info.json

[download](#)

```
ampliconsize: 2500
schemeversion: v1.0.0
schemename: artic-inrb-mpox
primer_bed_md5: 41dd2e95cbe2bda6be8556dd118c84fe
reference_fasta_md5: 7ca540c38981c844710721e5165d2b77
status: validated
citations:
authors: ARTIC network, INRB, Quick Lab
algorithmversion: primalscheme3
species: 10244
license: CC BY-SA 4.0
primerclass: primerschemes
infoschema: v2.1.0
articbedversion: v3.0
collections: QUICK-LAB, WHOLE-GENOME, ARTIC
links: { "protocols": [], "validation": [], "homepage": [], "vendors": [], "misc": [] }
refselect: { "KJ642613.1_masked": { "md5": "58062d2989d1d040b7fa06050a8b6040", "filename": "KJ642613.1_masked_refselect.fasta" } }
description: Clade I genomes were aligned to clade I reference and Clade II genomes were aligned to clade II reference with squirrel. The two clades were consensus aligned using MAFF
derivedfrom: null
contactinfo: null
```

<https://labs.primalscheme.com/detail/artic-inrb-mpox/2500/v1.0.0/?q=>

# Amplicon schemes dilutions

- Typically diluted to a working concentration of 10  $\mu$ M – **but not always!!!!**

Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design
ARTIC / INRB	400 (wastewater surveillance) 2500 (clinical surveillance)	147	<30	I and II (KJ642613.1)


When creating your primer pools check:

1. Primer dilution - by volume or concentration?
2. Pools – does the pooling parameters change? E.g Pan scheme vs Genotype specific PCR
3. What the stock concentration vs working concentration.
4. Clearly label your stock and working tubes!

# Amplicon Size

POOL 1

POOL 2

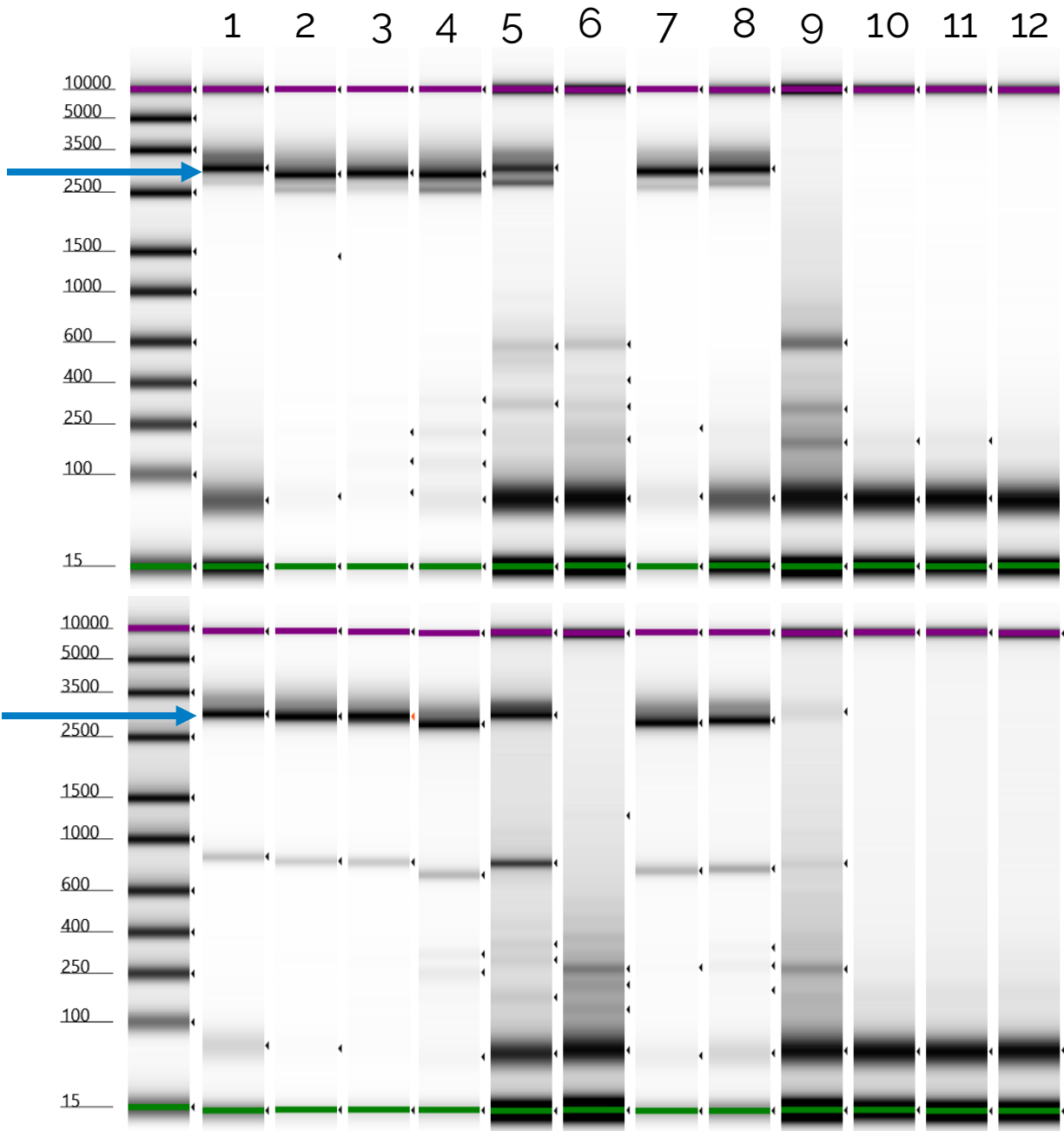


Sample	Mpox Ct
1	25.6
2	21.9
3	19
4	21
5	29.1
6	Not Det
7	24
8	26
9	35.3
10	Not Det
11	Not Det
12	Not Det

# Amplicon Concentration

POOL 1

POOL 2

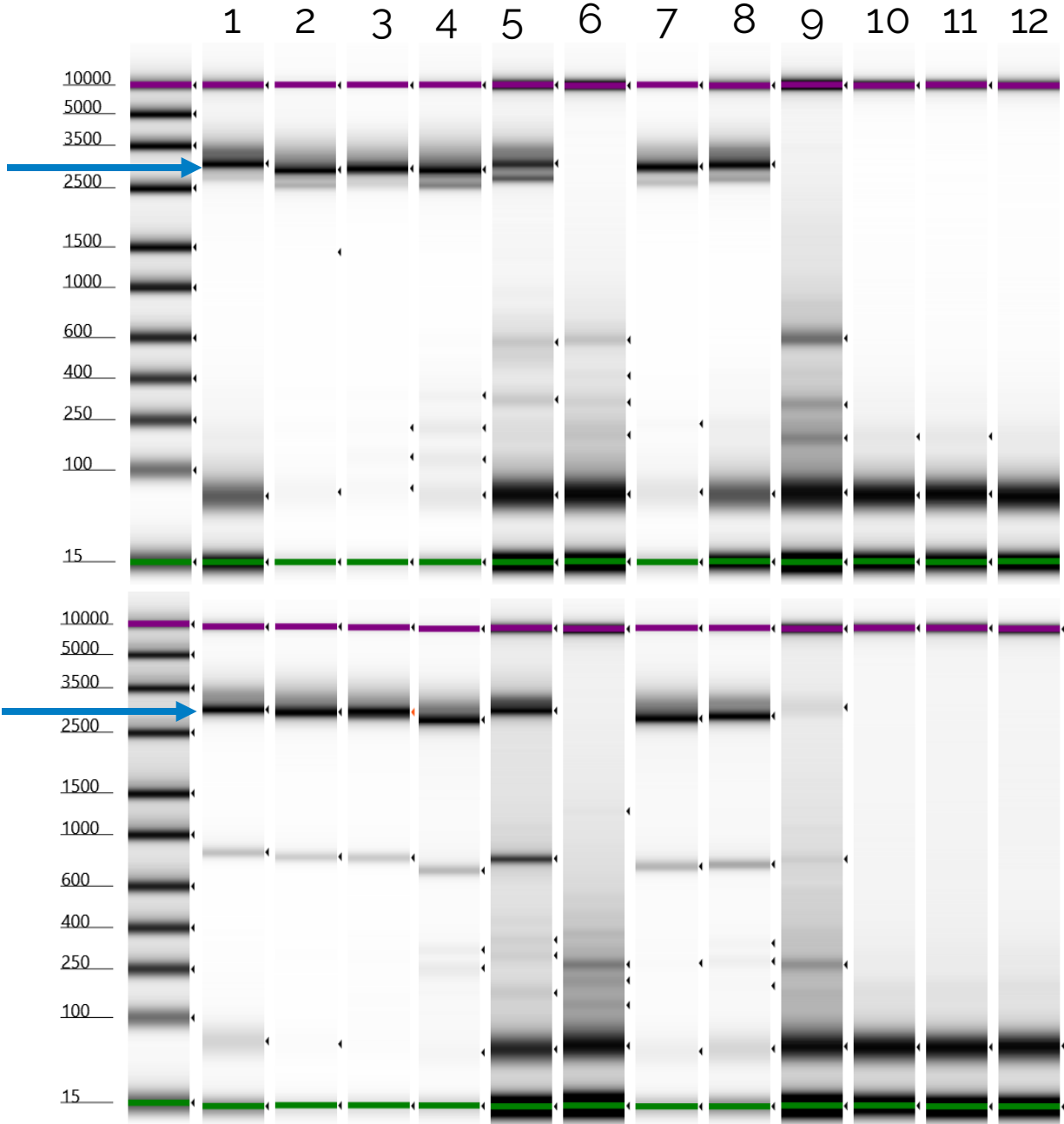

## QuBit Assay

Sample	Mpox Ct	POOL (ng/ul)	
		Pool 1	Pool 2
1	25.6	22.5	34.6
2	21.9	73.4	57.7
3	19	82.8	76.2
4	21	42.9	61.2
5	29.1	14.5	13.8
6	Not Det	13.8	10
7	24	42.2	49
8	26	25.3	31.2
9	35.3	13.8	14.3
10	Not Det	6.98	6.21
11	Not Det	9.08	50.6
12	Not Det	7.48	6.47

# Amplicon Summary

<https://www.nature.com/articles/nrmicro.2016.182#Sec6>

Method	Advantages	Disadvantages
PCR amplification sequencing	<ul style="list-style-type: none"><li>• Tried and trusted well-established methods and trained staff</li><li>• Highly specific; most sequencing reads will be pathogen-specific, which decreases sequencing costs</li><li>• Highly sensitive, with good coverage even at low pathogen load</li><li>• Relatively straightforward design and application of new primers for novel sequences</li></ul>	<ul style="list-style-type: none"><li>• Labour-intensive and difficult to scale for large genomes</li><li>• Iterating standard PCRs across large genomes requires high sample volume</li><li>• PCR reactions are subject to primer mismatch, particularly in poorly characterized or highly diverse pathogens, or those with novel variants</li><li>• Limited ability to sequence novel pathogens</li><li>• High number of PCR cycles may introduce amplification mutations</li><li>• Uneven amplification of different PCR amplicons may influence minor variant and haplotype reconstruction</li></ul>

# Wet Lab schedule DAY 1

MONDAY	ACTIVITY	PRESENTER
8:45 – 9:00	Registration	Lisa
9:00 – 9:15	Overview of the Doherty Institute/CPG/VIDRL/MDU	Lisa
9:15 – 9:45	Welcome and introductions	
9:45 – 10:00	Training overview	Jean
10:00 – 10:30	LAB: Pipetting exercise	Louise
10:30 – 11:00	Morning tea	
11:00 – 11:30	LECTURE: Introduction to MPXV genomics and VIDRL	Jean
11:30 – 12:30	LECTURE: Tiled Amplicon for Mpox	Jean
<b>12.30 – 13.30</b>	<b>Lunch</b>	
13.30 – 15.30	LAB: Tiled amplicon PCR	Louise
15.30 – 16:00	Afternoon tea	
16:00 – 16.30	LECTURE: Introduction to ONT sequencing viruses	Louise
16.30 – 17:00	Group discussion: Opportunity for Q&A and further discussion	Nicole