

## PRESENTER

### MONDAY

### ACTIVITY

8:45 – 9.00	Registration	Lisa
9.00 – 9.15	Overview of the Doherty Institute/CPG/VIDRL/MDU	Lisa
9.15 – 9:45	Welcome and Introductions	
9:45 –10.00	Training Overview	Jean
10.00– 10.30	<b>LAB: Pipetting exercise</b>	Louise
10.30 – 11.00	Morning tea	
11.00 – 11.30	<b>LECTURE:</b> Introduction to Mpox and MPXV genomics at VIDRL	Jean
11.30 – 12.30	<b>LECTURE:</b> Tiled amplicon Mpox viral theory	Jean
12.30 – 13.30	Lunch	
<b>13.30 – 15.30</b>	<b>LAB: TILED AMPLICON GENERATION</b>	<b>Louise</b>
15.30 – 16:00	Afternoon tea	
16:00 – 16.30	<b>LECTURE:</b> Introduction to ONT sequencing viruses	Louise
16.30 – 17:00	Group discussion: Opportunity for Q&A and further discussion	Nicole

# WET LAB DAY 1: TILED AMPLICON GENERATION

---

Dr. Louise Judd  
Centre for Pathogen Genomics- Innovation Hub  
Lead Laboratory Scientist and Genomics Trainer

# LAB REQUIREMENT S

# Avoiding Contamination: Lab Design

*Library preparation and Illumina sequencing by synthesis involve PCR amplification. It's important to avoid contamination to prevent the wrong template from being amplified which leads to inaccurate results.*

- Separate laboratory spaces for pre- and post- PCR steps with dedicated equipment and PPE in each
- At a minimum, perform pre-PCR in a sterile cabinet (PCR, laminar flow or biosafety)
- Work in a unidirectional manner
- We do not have the resources to run the training this way, but this is a very strong recommendation for your lab set up

# Avoiding Contamination: Pre- & Post- PCR

## Pre-PCR

### Master Mix Room

- For reagent preparation only
- Samples and nucleic acid should not be handled here

### Sample Preparation Room

- Samples and nucleic acid should be handled here
- Nucleic acid extractions performed here (kits stored here)

## Post-PCR

### Amplification Room

- Thermocyclers

### Post-Amplification Room

- Post-PCR analysis equipment
- Sequencers



# Avoiding Contamination: Lab Work

- Have a dedicated fridge/freezer for samples and a separate one for reagents --> or at least dedicated spaces
- Aseptic cleaning practises before and after work
  - 70-80% ethanol (EtoH), 5-15% Sodium Hypochlorite(NaOCl), DNAase Away
- Only use barrier tips to avoid pipette contamination
- Always use fresh tips for each pipetting step and for each different sample
- Use correct pipetting techniques
- Pulse-spin tubes and plates before opening to collect contents at bottom
- Open tubes and plates carefully, making sure to avoid any splashes
- Seal all plates and recap all tubes immediately after use

# Avoiding Contamination: NTC & PTC


*Remember to treat your controls in the same manner you treat your samples (i.e. if you add a reagent to your sample, add the same reagent to your controls). Try to follow the same control samples through your entire workflow.*

## NTC (Negative Template Control)

- Do not add any template DNA
- Can show:
  - reagent contamination
  - cross-contamination from samples

## PTC (Positive Template Control)

- Use a sample known to give a specific result
- Shows if our reagents and assays are working as expected

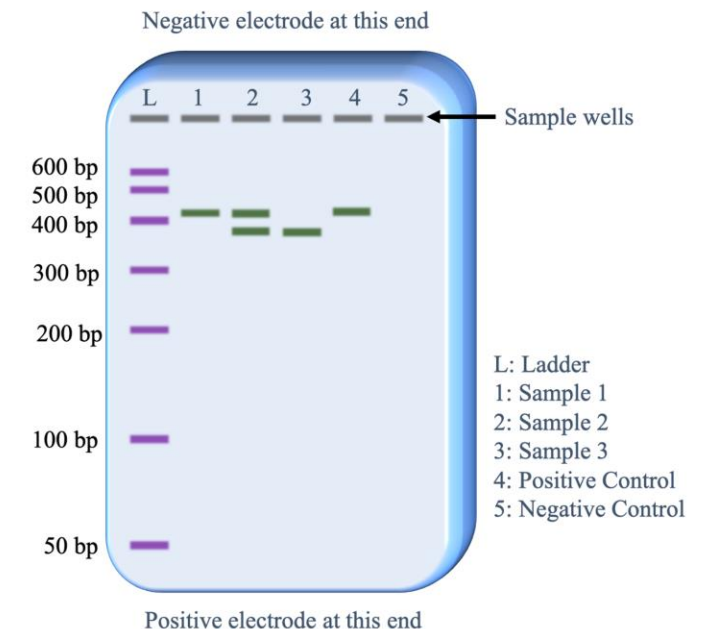

# Mpox SEQUENCING




# Mpox DNA requirements for Mpox sequencing

- DNA extracted from clinical swabs
- Quantitation of DNA is NOT required and is not informative
  - DNA from host and other microbes will be extracted alongside the Mpox DNA
  - The concentration of the extracted DNA does not correlate with the viral load in the sample
- Mpox is detected in the samples and viral load approximated with Mpox QPCR
- **Samples should have a Mpox QPCR cycle threshold (Ct) less than 30 to proceed with sequencing**
- Minimum volume of DNA extract required for sequencing 10 $\mu$ l

# Mpox DNA samples to be used for training

Sample number	Mpox Realtime Ct
1	25.6
2	21.9
3	Not detected
4	21
5	29.1
6	19
7	Not detected

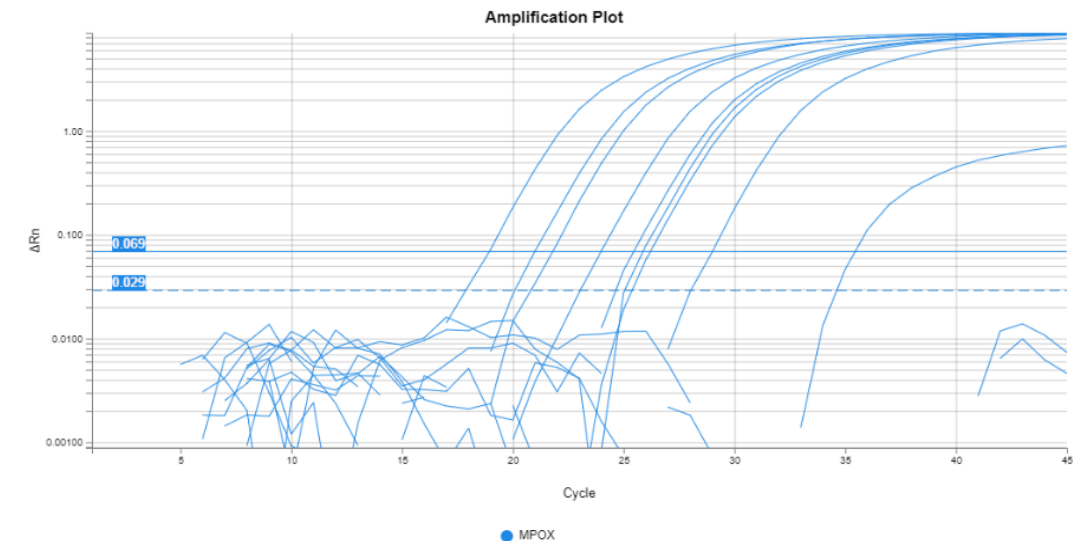

# PREPARING TILED AMPLICONS

# Mpox Panel for trainees

## a) Specimen collection - MPXV Panel for trainees.

PANEL OF CLINICAL SWABS RECEIVED AT VIDRL FOR MPXV TESTING

DNA EXTRACTED FROM SWAB

REAL TIME PCR SCREEN

# Extracting Mpox DNA

Kits we have tested

## Automated:

- Qiagen EZ2
- **MagMAX Viral/Pathogen Nucleic Acid Isolation Kit (MVP II)**
- QIAamp DNA QIAcube Kit
- TANbead OptiPure Viral Kit

## Manual:

- Qiagen QIAamp MinElute Virus Spin Kit

# Mpox Amplicon Scheme

Scheme	Amplicon size [bp]	# Primer pairs	Optimal Ct	Reference for primer design
ARTIC / INRB	400 (wastewater surveillance) <b>2500 (clinical surveillance)</b>	147	<30	I and II (KJ642613.1)

Pool 1 = 76 primers

Pool 2 = 71 primers



PCR 1 strip & PCR2 Strip  
(Pool 1 & Pool 2)

# Use of Strip tubes


## \*\*\*\*\*IMPORTANT REMINDER\*\*\*\*\*

- Always number the first of your tubes (the tube NOT the lid) of your strip with a number 1
- Strip tubes are symmetrical, it is visually not possible to know which is tube 1 and which is 8
- It is a very common error to flip the strip around during an experiment and get samples back to front
- Please make sure you always label the first tube is some way so you always know the correct orientation of your strip

# Use of strip tubes

## \*\*\*\*\*IMPORTANT REMINDER\*\*\*\*\*

Always number the first of your **tubes** (the tube NOT the lid) of your strip with a number 1

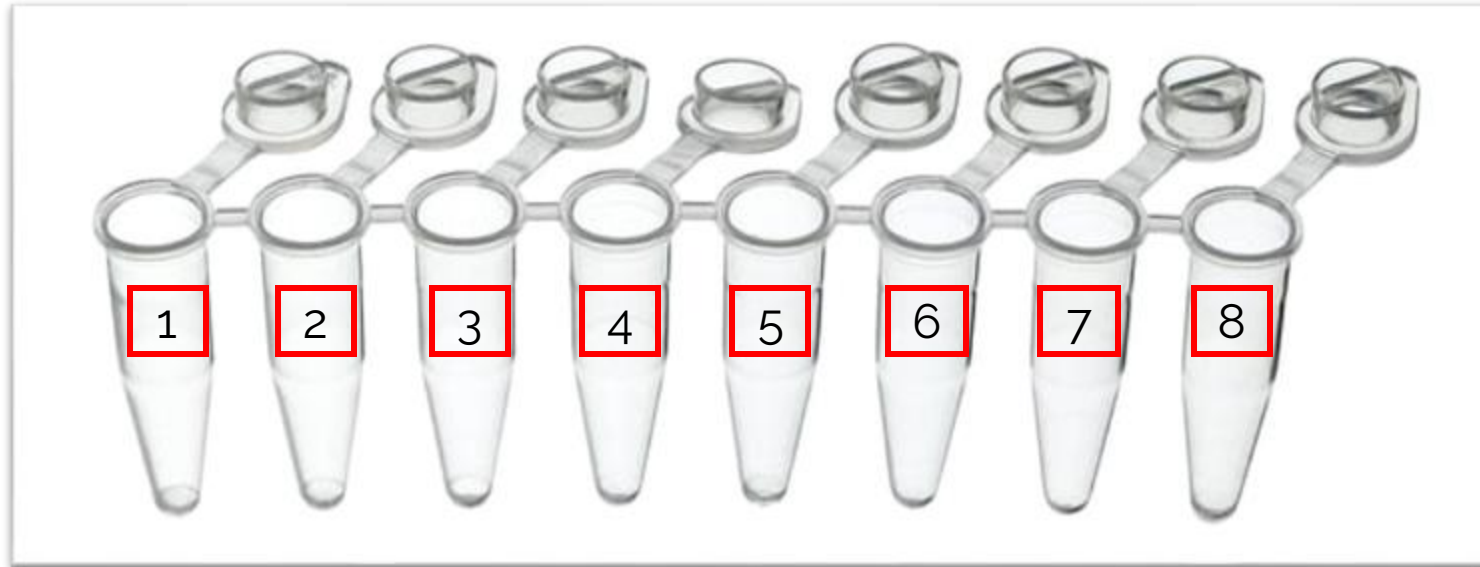



# Amplicon PCR Overview



## Reagents

Q5 Hot Start HF 2 x Master mix

*New England Biolabs M0494*


## Primer Pools 1 and 2

*INRB 2500 bp Amplicon Scheme*

*Stock 100  $\mu$ M- dilute 1:10 to generate working stock at 10  $\mu$ M*

Nuclease-Free Water (NFW)

2 x PCR strips (PCR1 and PCR2)

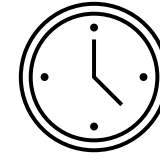

## Time

~ 4 hours

~ 30 minutes hands on

# Amplicon PCR- pool 1 mastermix

1. Label an Eppendorf tube "**PCR1 MM**" for your Amplicon Pool 1 Master Mix
2. Combine the following reagents in the PCR1 MM tube (amounts listed already account for overage)

Reagent	Single reaction	For 8 reactions (x9)
Nuclease free water	7.1 $\mu\text{l}$	63.9 $\mu\text{l}$
Mpox primer pool 1 (10 $\mu\text{M}$ )	2.9 $\mu\text{l}$	26.1 $\mu\text{l}$
Q5 Master mix (2x)	12.5 $\mu\text{l}$	112.5 $\mu\text{l}$
<b>Total</b>	<b>22.5 <math>\mu\text{l}</math></b>	<b>202.5 <math>\mu\text{l}</math></b>

3. Pipette to mix

# Amplicon PCR- pool 2 mastermix

1. Label an Eppendorf tube "**PCR2 MM**" for your Amplicon Pool 1 Master Mix
2. Combine the following reagents in the PCR1 MM tube (amounts listed already account for overage)

Reagent	Single reaction	For 8 reactions (x9)
Nuclease free water	7.3 $\mu\text{l}$	65.7 $\mu\text{l}$
Mpox primer pool 2 (10 $\mu\text{M}$ )	2.7 $\mu\text{l}$	24.3 $\mu\text{l}$
Q5 Master mix (2x)	12.5 $\mu\text{l}$	112.5 $\mu\text{l}$
<b>Total</b>	<b>22.5 <math>\mu\text{l}</math></b>	<b>202.5 <math>\mu\text{l}</math></b>

3. Pipette to mix
4. Clearly label 2 strip tubes as "PCR 1" and "PCR 2" with your name and initials

# Amplicon PCR – pool 1 PCR

1. Transfer **22.5 µL of PCR1 MM** to each well in the **PCR1** strip tube
2. Add 2.5 µL of DNA sample 1 to the well 1 of PCR1strip tube
3. Mix gently by pipetting up and down
4. Change pipette tips between samples
5. Add 2.5 µL of DNA sample 2 to the well 2 of PCR1strip tube
6. Repeat this with remaining DNA samples 3-8 into appropriate well
7. Ensure each tube is capped
8. Spin down briefly

# Amplicon PCR –pool 2 PCR

1. Transfer **22.5 µL of PCR2 MM** to each well in the **PCR2** strip tube
2. Add 2.5 µL of DNA sample 1 to the well 1 of PCR2 strip tube
3. Mix gently by pipetting up and down
4. Change pipette tips between samples
5. Add 2.5 µL of DNA sample 2 to the well 2 of PCR2 strip tube
6. Repeat this with remaining DNA samples 3-8 into appropriate well
7. Ensure each tube is capped
8. Spin down briefly

# Amplicon PCR – running PCR

1. Run the Mpox PCR program on the thermocycler

## **Mpox PCR Program**

Preheat Lid: 105°C

Reaction Volume: 25 µL

.....

98°C for 30 secs minutes

35 Cycles of:

98°C for 15 seconds

65°C for 5 minutes

Hold at 10°C

# Amplicon PCR – running PCR

1. Run the Mpox PCR program on the thermocycler

## Mpox PCR Program

Preheat Lid: 105°C

Reaction Volume: 25 µL

.....

98°C for 30 secs minutes

35 Cycles of:

98°C for 15 seconds

65°C for 5 minutes

Hold at 10°C

Run PCR Program overnight


PCR 1 strip & PCR2 Strip  
(Pool 1 & Pool 2)

## Day 2 Qubit QC of Amplicons



# THANK YOU

[centre-pathogen@unimelb.edu.au](mailto:centre-pathogen@unimelb.edu.au)

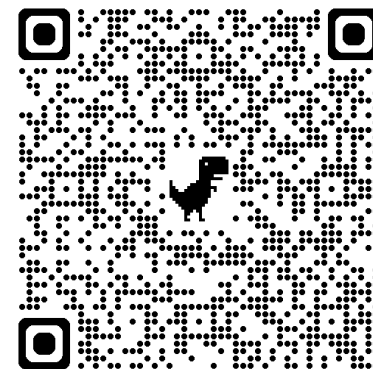



## PRESENTER

### MONDAY

### ACTIVITY

8:45 – 9.00	Registration	Lisa
9.00 – 9.15	Overview of the Doherty Institute/CPG/VIDRL/MDU	Lisa
9.15 – 9:45	Welcome and Introductions	
9:45 –10.00	Training Overview	Jean
10.00– 10.30	<b>LAB:</b> Pipetting exercise	Louise
10.30 – 11.00	Morning tea	
11.00 – 11.30	<b>LECTURE:</b> Introduction to Mpox and MPXV genomics at VIDRL	Jean
11.30 – 12.30	<b>LECTURE:</b> Tiled amplicon Mpox Viral Theory	Jean
12.30 – 13.30	Lunch	
13.30 – 15.30	<b>LAB:</b> Tiled Amplicon Generation	Louise
<b>15.30 – 16:00</b>	<b>AFTERNOON TEA</b>	
16:00 – 16.30	<b>LECTURE:</b> Introduction to ONT sequencing viruses	Louise
16.30 – 17:00	Group discussion: Opportunity for Q&A and further discussion	Nicole